Richtanalyse in Masse-%

Elemente		%	
Pd	57,5	Ag	31,5
In	1,9	Ru	
		Sn	9,0

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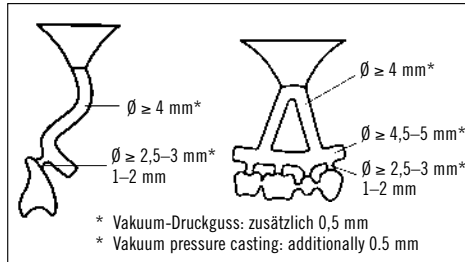
**Instructions for use****Modelling:**

- Minimum metal thickness (after grinding): for ceramic veneering 0.4 mm, for acrylic veneering with retention pearls 0.3 mm.
 - Connecting parts between the pontics should be as thick and high as possible (at least 3.5 mm high and 2.5 mm wide).
 - Spruing of single crowns: provide casting reservoir.
- Investing:** Use phosphate-bonded crown and bridge investment materials

(e.g. Bellavest®, BellaStar).

Casting and finishing: General: Do not overheat alloy. Use only clean crucibles, one crucible per alloy. Recommendation: to enable an exact identification of each case cast new metal only.

- In case of re-casting: only re-cast identical alloys. Blast old material. Add at least 50 % of new material.
- Use ceramic crucible and sprinkle some Auromelt HF melting powder on the casting ingots.
- Continue to heat after complete melting of the casting pieces:

Flame melting:	10–15 seconds
HF induction heating:	approx. 10 seconds
Resistance heating:	2–3 minutes.

 Nautilus T: 6213 to 6216 (ceramic crucible/glass carbon cylinder)
 Nautilus CC plus: 172 (ceramic crucible)
 165 (glass carbon cylinder)

Use fine carbide, ceramically bonded stones or BEGO sintered diamond milling tools for finishing.

Ceramic: Only use veneering ceramics which are in accordance with ISO 9693 (e.g. VM13, Inline) with firing temperatures of up to approx. 980 °C and which are resistant to green discoloration by means of silver oxides! Always follow the ceramic manufacturer's instructions!

- Always blast the surface to be veneered (Korox® 110, 2 – max. 3 bar) and clean the frame thoroughly (steam clean or boil in aqua dest.).
- Long-term cooling recommended (cooling phase up to approx. 600 °C).

Oxide firing:

- Maintain oxide firing at 960 °C without vacuum for 10 minutes.
- The oxide shall be blasted again prior to application of the ceramics.

Soldering:

- Support object in a soldering block of Bellatherm®. Prepare a gap of max. 0.2 mm with parallel walls.
- Soldering before firing with the flame (1125 °C): BegoStar®-Solder (REF 61081) and Minoxid or Fluxsol.
- Soldering after firing in furnace (810 °C): BEGO-Gold-Solder I (REF 61017) and Minoxid. Allow to cool normally.
- Acid-treat residual flux in Aurocid (60 °C, approx. 1 minute). Clean thoroughly (steam clean or boil in aqua dest.).

Laser welding: Filler material: BegoPal® 300-wire diameter 0.35 mm (REF 61165) or accurately fitting, custom-ground casting pieces.

Prescription only / For professional use only

Secondary effects: Such as allergies to contents of the alloy or electrochemically based reactions may very rarely occur.

Reciprocal actions: In case of occlusal or approximal contact of different alloys electrochemically based reactions may very rarely occur.

Reactions: In case of known incompatibilities and allergies to contents of the alloy.

Warranty: Whether given verbally, in writing or by practical instructions, our recommendations for use are based upon our own experience and trials and can only be considered as standard values. Our products are subject to a constant further development. Therefore alterations in construction and composition are reserved.